



APP



SHUNMEI MEDICAL CO., LTD.

Add: R301 and R401, Building B, Shenchanggang Investment Co., Ltd.,
No. 8, Jinlong First Road, Baolong Community, Baolong Street,
Longgang District, Shenzhen, 518116 Guangdong, P.R. China

Branch Office: Yifa 3rd road, Yifa Industrial Zone,
Pingtan Town, Huiyang District, Huizhou,
516259, Guangdong, P.R. China

Mobile/Whatsapp: +86 18344359973 +86 18344369007

E-mail: nora@shunmed.com nora1@shunmed.com

Website: www.shunmed.com




Shunmei

Search



Double Lumen Microcatheter

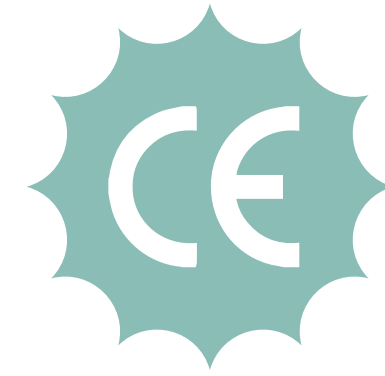


Precision
Twice
Reach
Further

Medical technology creates future

Shunmei®

Adantages:

- Double Lumen Catheter can accommodate two guidewires simultaneously while ensuring excellent controllability, support, and crossability.
- It can prevent the guide wires from entangling with each other.

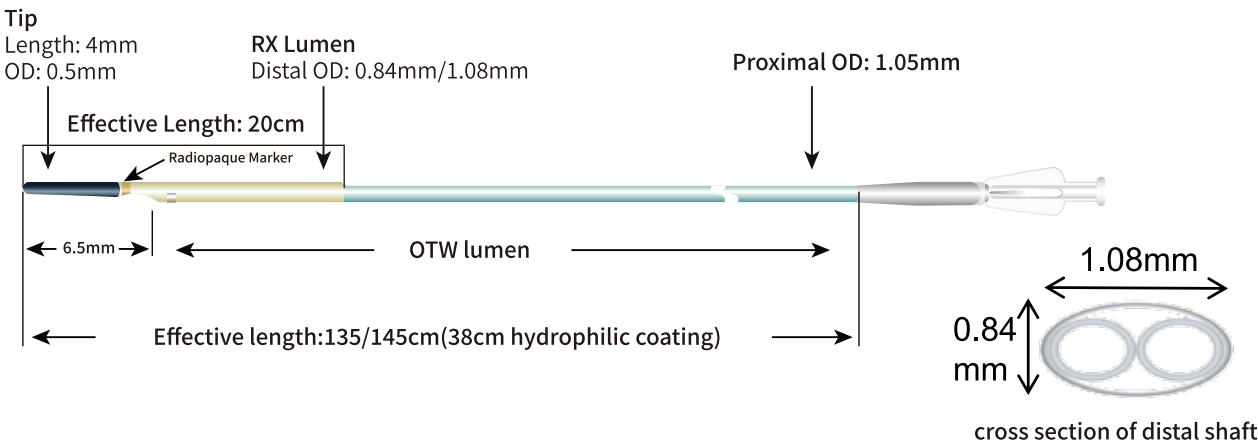

Product Specification

Code	Model	Maximum outer diameter of proximal catheter	Maximum outer diameter of distal cathete	Effective Length	Inner diameter of over the wire lumen(OTW)	Inner diameter of rapid exchange lumen(RX) / tip diameter
639950	DAMC-145	1.05mm	1.08mm	145cm	0.43mm	0.41mm
639951	DAMC-135	1.05mm	1.08mm	135cm	0.43mm	0.41mm

Product Features

6.5mm Tip To OTW Port

Provide better support and control for the OTW guide wire

Dual Radiopaque Rings

Precise Marking of Two Guidewire Exits

Excellent Crossability

The double-tapered design at the tip of the RX lumen offers excellent crossability and trackability

38cm Hydrophilic Coating

Exceptional Lubricity